


# LSRP-F-3

3 Cubic Ft. Undercounter Laboratory Freezer

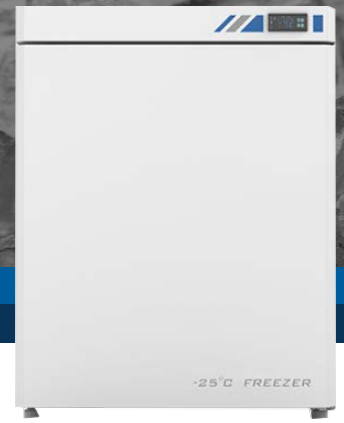

**Temperature Range**  
-10°C to -25°C

**Exterior W x D x H**  
27.5" x 26.3" x 33.6"

**Capacity**  
3 cu.ft.

**Weight (lb)**  
120 lbs

## Undercounter Laboratory Freezer

Undercounter laboratory freezer suitable for the storage of samples and reagents. Suitable for use in hospitals, clinics, and laboratories.

### Multiple Alarms Ensure Safety of Storage

The freezer includes audible & visual alarm functions, such as a high temperature/low temperature alarms, door ajar alarm, and sensor failure alarm to guarantee security of samples in a comprehensive way.

### High-precision Temperature Control

The high-precision computerized temperature control system ensures an adjustable temperature within a range from -10 to -25°C inside the cabinet, and the high-definition digital temperature display has a display accuracy of 0.1 °C.

### Insulation

The polyurethane insulated door delivers great thermal insulation to ensure a consistent internal temperature.

- Capacity 3 cu.ft.
- Internal Dims(w x d x h):  
19.6" x 19" x 21"
- External Dims(w x d x h):  
27.5" x 26.3" x 33.6"
- Weight: 120 lbs
- Temperature Range: -10C ~ -25C
- Ambient Temp: 16C - 32C
- Controller: Microprocessor
- Digital Display
- Manual Defrost
- Refrigerant: R600a
- External Material: Powder Coated
- Inner Material: Aluminum spray coated plates
- 3 Shelves
- Door lock with key
- Access Port: 1 Ø 25mm
- Casters/Levelers: 2 / 2
- Alarms: High/Lo Temp, Sensor Failure, Door ajar
- Power Supply: 115V/60Hz
- Current: 1.26 A



**This unit comes with a standard 2 year warranty.**